


## Waterbath

### WNB 22

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at [myAtmoSAFE@memmert.com](mailto:myAtmoSAFE@memmert.com).

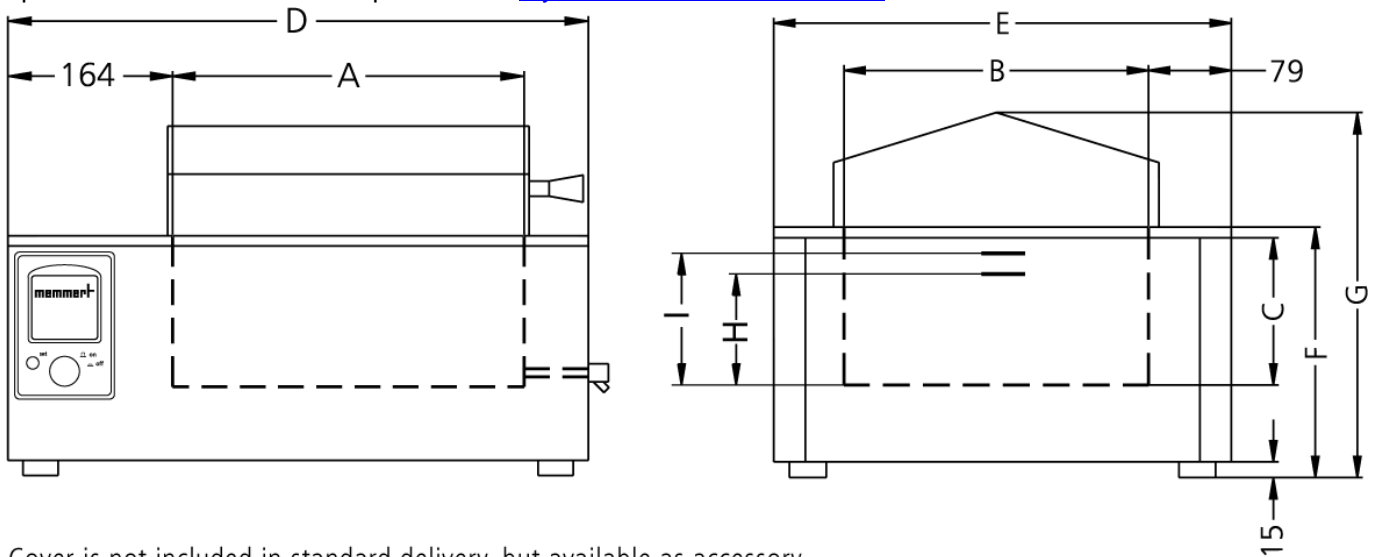

Cover is not included in standard delivery, but available as accessory.

## Temperature

<b>Setting temperature range</b>	+10 to +95 °C and boiling stage
<b>Working temperature range in °C</b>	min. 5 above ambient up to +95 °C with additional boiling mode

## Control of standard components

<b>Temperature sensor</b>	1 Pt100 sensor class A in 4-wire-circuit
<b>Timer</b>	integrated digital timer from 1 min. to 99,99 hours for:ON continuous operationWAIT (delayed on for continuous andlimited timed operation)HOLD

## Safety

<b>Temperature control</b>	mechanical temperature limiter TB protectionclass 1 switching the heating off at approx. 30°Cabove max. temperature of the bath
<b>Temperature control</b>	in case of overtemperature due to failure,the heating is switched off at approx. 10°Cabove the set temperature (fixed value)
<b>Autodiagnostic system</b>	microprocessor PID-temperature controllerwith integrated autodiagnostic system withfault indicator

## Heating concept

<b>Heating Baths</b>	corrosion-proof large-area heating onthree sides
----------------------	--

## Standard equipment

<b>Controller</b>	digital display (LED) of set and actualtemperature (0,1°C resolution) and of (remaining)programme time
-------------------	--

## Stainless steel interior

<b>Dimensions W x H x D in mm</b>	$w_{(A)} \times h_{(C)} \times d_{(B)}$ : 350 x 220 x 290 mm
<b>Interior</b>	easy-to-clean interior, made of stainlesssteel, reinforced by deep drawn ribbing,material no. 1.4301 (ASTM 304), laser-welded
<b>Volume</b>	22 l
<b>Liquid level min.</b>	(H) 177 mm
<b>Liquid level max.</b>	(l) 200 mm

## Textured stainless steel casing

<b>Dimensions</b>	$w_{(D)} \times h_{(G)} \times d_{(E)}$ : 578 x 405 x 436 mm
-------------------	--

### Electrical data

**Voltage** 230 V, 50/60 Hz

**Electrical load** approx. 2000 W (during heating)

### Ambient conditions

**Installation** The vent openings in the left and back side must remain unobstructed. Minimum wall spacing on all sides is 80 mm. The minimum spacing from the top of the bath to the next ceiling is 750 mm.

**Ambient temperature** +5 °C to +40 °C

**Humidity rh** max. 80 %, non-condensing

**Overvoltage category** II

**Pollution degree** 2

### Packing/shipping data

**Transport information** The appliances must be transported upright

**Customs tariff number** 8419 8998

**Country of origin** Federal Republic of Germany

**WEEE-Reg.-No.** DE 66812464

**Dimensions approx incl. carton** w x h x d: 670 x 530 x 460 mm

**Net weight** approx. 16 kg

**Gross weight carton** approx. 23 kg

**Standard units are safety-approved and bear the test marks**

