



# B6e

## THREE FUNCTIONS ELECTRIC BED

Registration code: SKD-C







**SAIKANG** / Medical

Healthier  
your life

# B6e

## FEATURES AND HIGHLIGHTS

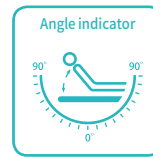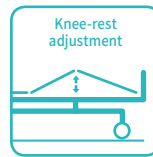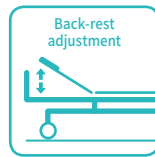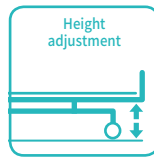


### Technical parameters

External size (LxW):	2170x1060 mm
Mattress platform:	1925x900 mm
Safe working load:	250 kg
Height adjustment:	435-745 mm
Back-rest adjustment:	0-75° (±10°)
Knee-rest adjustment:	0-35° (±10°)
Angle indicator:	0-90°
Voltage and frequency:	100-240 V · 50/60 Hz
Plug as required (hospital grade)	




### Technical configuration

<input checked="" type="checkbox"/> Linear Motor	3 pcs
<input checked="" type="checkbox"/> Hand Controller	1 pc
<input checked="" type="checkbox"/> 5" Hardshell Castors	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail Extra Height	1 set
<input checked="" type="checkbox"/> ABS Bed Ends with Safe Lock	1 set
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	2 pcs

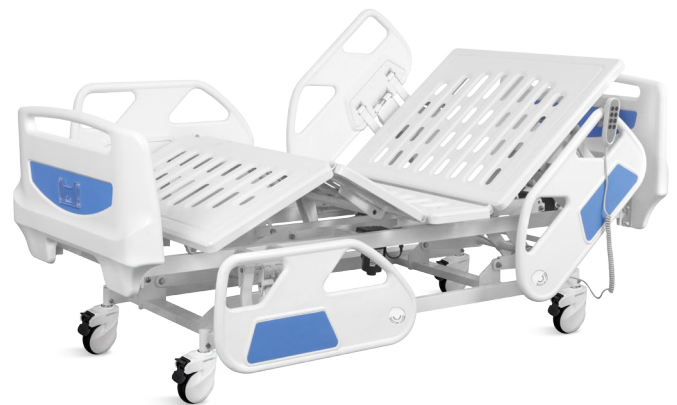


### Optional

<input type="checkbox"/> Manual CPR	
<input type="checkbox"/> Electric CPR	
<input type="checkbox"/> Back-up Battery	
<input type="checkbox"/> X-ray Bed Surface	
<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan Traction Frame	
<input type="checkbox"/> Directional Wheel	

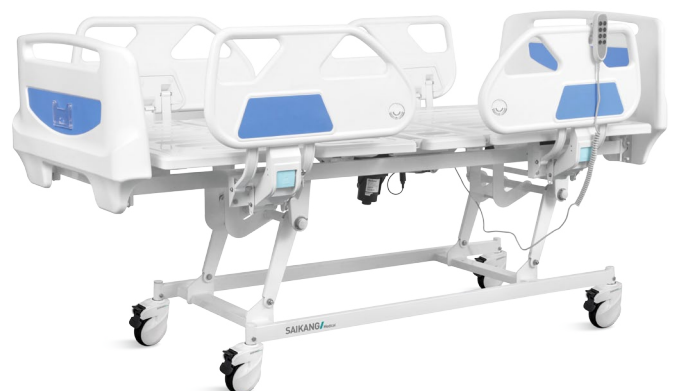

**Design**

E-shaped large guardrail, with angled display and built-in controller, meets safety standards, providing comforting feeling for patients.

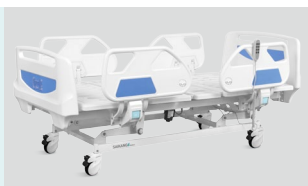**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**


Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Angle Indicator**

Patented angle indicator.

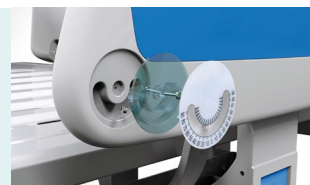**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Side Rail**

Full ABS material, tuck away side rail with built-in controller. From top of the guardrail to bed surface is about 345mm, more height provides better protection for patient. Bottom can endure 100,000 times press, buffer, built-in angle indicator.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000 times press. Head and foot board can be reversible.

**Wheels**

TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**IV pole**

Stainless Steel Column with plastic hook.

**Drainage Hooks**

Mobile plastic drainage hooks.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue