



Grace Care series

# D3w

## THREE FUNCTIONS MANUAL BED

Registration code: SK-CD1




**SAIKANG** Medical

Healthier  
your life

# D3w

## FEATURES AND HIGHLIGHTS

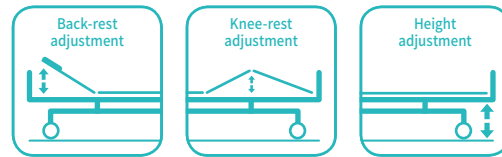


### Technical parameters

▪ External size (LxWxH)	2245x1000x435-745 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)



### Technical configuration

<input checked="" type="checkbox"/> Manual Cranks	3pcs
<input checked="" type="checkbox"/> 5" Hardshell Castors	4pcs
<input checked="" type="checkbox"/> Foldable Side Rail with ProErgo-Handle™	1set
<input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> IV Pole	1set
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs

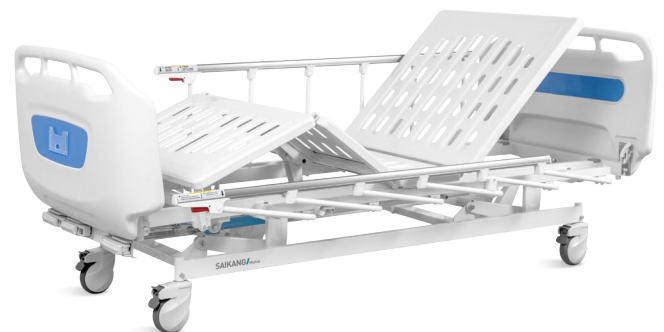

### Optional

<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan Traction Frame	
<input type="checkbox"/> Directional Wheel	



**Design**

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa  
Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11  
Process epoxy painting, ASTM testing  
anti-bacterial, paint thickness 0.12  
mm, brightness 60°, paint can resist  
50 kg impact.

**Welding | Chassis**


Panasonic robotic ensure, 360° full  
smooth welding Solid bed base  
made of 5.0 mm stamped rotating  
plate.

**Hand Control**

Retractable hand crank, 100,000  
times remaining intact  
Extreme protection device, sounds  
when tipped to the max position,  
Strong pallets to increase the bed  
height adjustment (435-745mm).

**Platform**

Stamping forming platform, smooth  
surface no welding scar, bearing  
250KG load. Option resin X-ray  
translucent platform.

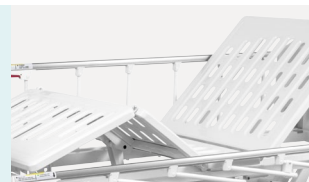**Guardrail**

L-shaped foldable guardrail (Patent)  
with ProErgo-Handle™, patented parts.  
with red and green mark to indicated  
if the siderail is locked or not. Bended  
tube, solid installing, 0 shaking, anti-  
pinch, protect patient's safety.

**Head & Foot Board**

Head and foot board can be reversable,  
extra height design, size: 1000\*563mm;  
ABS Grace Bed Ends with Safe Lock,  
automatic lock without any manual  
operation on lock, on/off removable  
mark on hook lock.

**Castors**

5" hardshell casters  
TPR tire no worn out after running  
30KM, save anti-winding hard shell,  
united forming without bolts.  
Pass Dynamic test: Bearing 120kg run  
30KM and pass obstacles 500 time.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic  
hook.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue